Instructions Manual Heating Pad

Model: B-HT211

BEIGE | BLUE

Thank you for choosing the Heating Pad. Our goal is to provide the best products to enhance your life.

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WARNING!

RISK OF ELECTRICAL SHOCK

The product is designed specifically for humans.

WARNING: This is an electrical product. Please read all directions before use. Failure to do so may result in electrical shock or fire causing serious bodily injury.

- Read all instructions before using this product.

- FOR INDOOR USE ONLY

- DO NOT use if any visible damage is apparent (cracks, tears or exposed wires and heat element).

- Must instantly turn off the Heating pad and unplug when overheat or charred flavor in use.

- Serious heart patients, those that depend on a pace maker and pregnant women should check with their doctor before use.

- DO NOT use for unconscious and sensitive to heat people.

- DO NOT make any modifications to switch button. wires and heat element.

- DO NOT use in sleeping; Do not put the heating pad under in mattress or quilt when use.

- DO NOT use cord as handle.

- DO NOT puncture or use any sharp object to secure the Heating pad.

- DO NOT use another electronic heating device in conjunction with your heating pad.

- DO NOT fold the heating pad when use.

- The maximum weight which can sit and lay on the heating pad is 100 kg. Do not sit and sleep on when the heating pad is folded.

- DO NOT Heap and Cover anything on the heating pad when use.

- DO NOT USE for Children.

- Unplug when not in use.

- DO NOT USE WET.

- DO NOT dry clean this heating pad. Dry cleaning solvents may damage the wires.

- DO NOT use bleach or cleaning fluids.

- DO NOT use a wringer.

- DO NOT iron your heating pad.


- DO NOT use mothballs or sprays of any kind.

- DO NOT reconnect the heated product or use it until it is COMPLETELY dry.

- DO NOT use analgesic cream in combination with your heating pad.

Function:

- Offers excellent coverage for large muscle groups

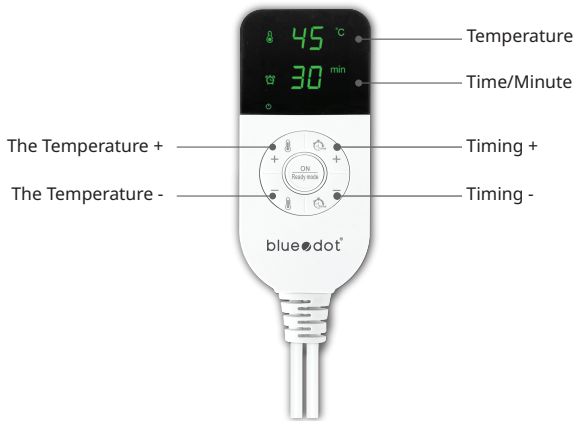- Temperature Settings: LED controller adjusts heat easily with temperature settings to accommodate your comfort level (30°C to 70°C).
- TIMER: Settings range from 10 to 90 min, conserves energy and prevents burns.
- 200 cm Extra-Long Cord: Makes it easier to use from the comfort of your bed or favorite chair.
- Machine Washable: Helps maintain the beauty and softness for years to come.
- Clean: Turn off the Heating Pad, Unplug heat pad from the wall outlet. Remove the controller and Power unit. heat pad can be machine washed and must be dry before reassembling the unit. Do not use cleaners or chemicals.
- Storage: Turn off the heat pad, dial out the heat pad, Store the pad in a dry, safe, clean and cool place after clean up.
- Full wire constant temperature system: NTC+PTC, with case temperature protection, and thyristor breakdown protection, full wire constant temperature protection, thermal fuse and current safety protection, 102°C thermal fuse and 1A current safety protection. 70°C power off function, can be resumed to work when temperature going down below 70°C.

SAFETY:**CAUTION**

The heating pad is fitted with overheating protection SAFETY SYSTEM.

The electronic sensor technology prevents the entire surface of the heating pad from overheating by automatically switching off in the event of a fault. If the safety system switches the heating pad off automatically as a result of a fault, the temperature display on the control will no longer be illuminated when the heating pad is switched on.

Please note that for safety reasons, the heating pad must not be used after a fault and must be returned to the service center.

PRE-OPERATING INSTRUCTIONS

When removing the unit from the packing case, check for any shipping damage.


Please contact your shipper and dealer immediately if the unit bears any sign of damage from mishandling.

Contents:

1 x Heating pad

1 x controller box with power cord

1 x Instructions Manual

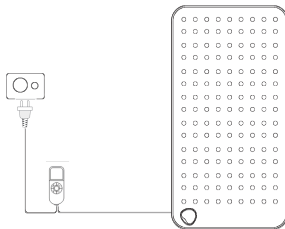**Limitations of Liability:**

1. Terms of usage is unconditional upon your acceptance, without changes to the terms, conditions or notice listed herein.
2. Proper use including reviewing all safety precautions.
3. No prohibited or unlawful use. The product is intended for humans use only.
4. Our reserves the right to change the terms and conditions under which this product is offered.
5. Buyer assumes all risks and liability from the use of this product. Our is not liable for any damages arising from use or misuse of this product.

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Product Owner: Lifebox Asia Company Limited